


Hill-Rom® Procedural Stretcher



NEW

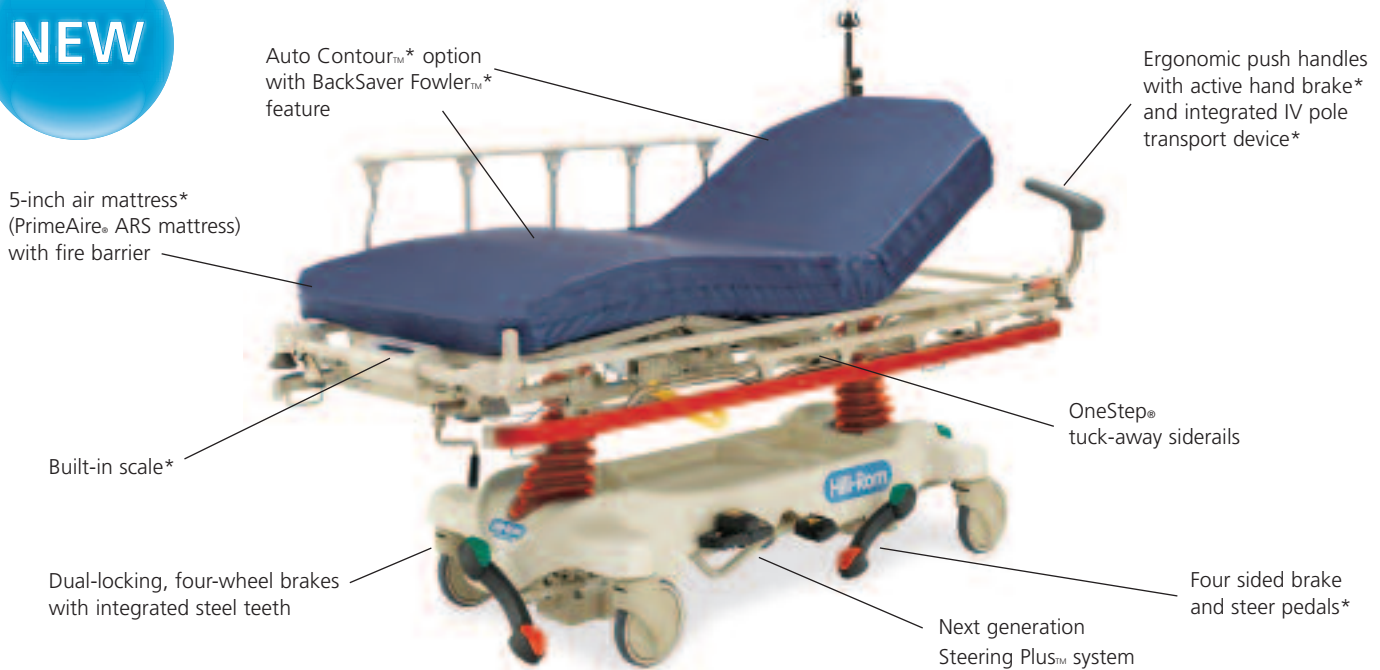

*Optional features

At a glance

New Hill-Rom stretcher key features:

- **700 lbs. weight capacity**
- **Built-in scale** shows kilograms and pounds simultaneously
- **Low height (20.7")** for safer stretcher ingress and egress – best in class
- **5-inch air mattress** (PrimeAire ARS® mattress) for comfortable patients – the first non-powered air mattress for a stretcher
- **BackSaver Fowler™ feature** for easier articulation of the head section
- **Auto Contour™ option** helps prevent patients from slipping down the mattress
- **Active hand brake** – the world's first stretcher hand brake to help you negotiate ramps and corners without using your body weight and for safe and easy stops
- **Integrated IV pole transport device** helps to maneuver stretcher and a portable IV pole with ease
- **Next generation Steering Plus™ system** for crisp cornering and control
- **Dual-locking, four-wheel brakes** for a strong, stable platform
- **Power washable**
- **Three year parts warranty** and life-time frame weld warranty

New Hill-Rom Procedural stretchers

Hill-Rom's rugged procedural stretcher (P8000) is the ideal choice for emergency, OR, PACU or outpatient settings. This stretcher offers outstanding mobility, comfort and quality.

Mobility. Why struggle to steer and stop? An optional Active Hand Brake reduces strain of going down ramps and sudden stops. Our integrated IV pole transport device helps you to maneuver stretcher and portable IV pole stands with ease. The next generation Steering Plus™ system allows crisp cornering and control.

Comfort. Reduce patient-handling risk and enhance your and your patient's comfort. The BackSaver Fowler™ feature provides you with easier articulation of the head section. Combine it with the Auto Contour™ feature to minimize repositioning of patients. Our new 5-inch air mattress (PrimeAire® ARS mattress) with fire barrier is the industry's first non-powered air surface for stretchers — helping to make patients more comfortable.

Quality. Our reinforced frames are made of heavy-gauge, automotive-grade steel for long-lasting quality – and a 700 lbs. weight capacity to handle higher patient weights. Nickel-plated push handles, IV poles, siderails and other components offer superior rust resistance, even in harsh environments.

Hill-Rom®

Hill-Rom® Procedural Stretcher



| Features | Standard Features | Optional Features |
|---|-------------------|-------------------|
| 700 lbs. weight capacity | • | |
| Built-in, dual-color display scale | | • |
| 5-inch PrimeAire® Air Reflex™ System (ARS) mattress with fire barrier | | • |
| BackSaver Fowler™ (includes manual knee gatch) | | • |
| Auto Contour™ option (includes "BackSaver Fowler") | | • |
| Active hand brake | | • |
| Integrated IV pole transport device | | • |
| Brake/Steer pedals – Head and foot | • | |
| Brake/Steer pedals – All four sides | | • |
| Next generation Steering Plus™ system | • | |
| Ergonomic push handles | | • |
| Permanent IV pole – 3 stage | | • |
| 26-inch width | • | |
| 30-inch width | | • |
| Manual knee gatch | | • |
| Footend hi/lo hydraulic pedals | | • |
| Siderail system | Tuck-away | |
| Radiolucent deck | | • |
| 3-inch high density foam mattress with fire barrier | • | |
| 4-inch high density foam mattress with fire barrier | | • |
| 5-inch Comfortline® multi-density foam mattress with fire barrier | | • |
| Bumper Colors – Choice of neutral, blue, red, teal or purple | | Choice |

- Accessories**
- Stretcher pillow
 - IV transporter
 - Patient tray
 - Utility tray
 - Restraint straps
 - Oxygen tank holder
 - Liquid oxygen tank holder
 - Infusion support system
 - Siderail covers
 - Armboard
 - Detachable footboard
 - Mattress pad
 - Lateral X-ray cassette holder
 - Chart holder
 - Convertible footboard
 - Accessory IV pole
 - Ankle stirrups
 - Transport straps
 - Foot extender
 - Paper roll dispenser

- Technical Data**
- Maximum weight limit:700 lbs. (317.5 kg)
- Low position –
 floor to top of surface:20.7" (52.6 cm),
 with scale: 20.5" (52.1 cm)
- High position –
 floor to top of surface:34.25" (87 cm)
 with scale:34" (86.4 cm)
- Overall length:83" (210.8 cm)
- Overall width (siderails up):32" (81.3 cm)
or 36" (91.4 cm)
- Overall width (siderails stored):
26.5" (67.3 cm) or 30.25" (76.8 cm)
- Siderail length:47" (119.4 cm)
- Siderail height above
 mattress/surface:14" (35.6 cm)
- Mattress size:
26" (66 cm) x 75" (190.5 cm)
or 30" (76.2 cm) x 75" (190.5 cm)
- Maximum head elevation:90°
- Maximum Trend/Reverse Trend:18°
- Floor to base clearance:3.5" (8.9 cm)
- Caster size:8" (20.3 cm)



Auto Contour™ option with BackSaver Fowler™ feature



Built-in Scale



5-inch Air Mattress (PrimeAire® ARS mattress)



Integrated IV Pole Transport Device



Active Hand Brake

Hill-Rom reserves the right to make changes without notice in design, specifications and models. The only warranty Hill-Rom makes is the express written warranty extended on the sale or rental of its products. © 2006 Hill-Rom Services, Inc. ALL RIGHTS RESERVED. 145647 12/21/06