

Hill-Rom®

Procedural Stretcher



At a glance

Hill-Rom stretcher key features:

- **700 lb weight capacity**
- **Low height** (20.7") for safer stretcher ingress and egress – best in class
- **Therapeutic Support Systems:** Foam, TempurPedic®, AccuMax Quantum*
- **BackSaver Fowler™ mechanism** provides for easier articulation of the head section*
- **Auto Contour™ feature** helps prevent patients from slipping down the mattress* (includes BackSaver Fowler mechanism)
- **Active hand brake** reduces physical exertion required when going down ramps, around corners and stopping*
- **Integrated dual display scale** (lbs/kgs simultaneously) memorizes last weight*
- **Ergonomic push handles***
- **Integrated IV pole transport holder** helps to maneuver stretcher and a portable IV pole with ease*
- **Integrated O₂ holder and tray***
- **Next generation Steering Plus™ system**
- **Redundant-locking, four-wheel brakes**
- **Power washable**
- **Three year parts warranty** and life-time frame weld warranty

* Optional feature

Hill-Rom Procedural Stretcher

Hill-Rom's rugged procedural stretcher (P8000) is the ideal choice for Emergency, OR, PACU or ASC settings. This stretcher offers unsurpassed mobility, comfort and quality.

Mobility. Hill-Rom's innovative Steering Plus™ system allows crisp cornering and control. In addition, the system self adjusts when passing over uneven floor surfaces, preventing vibrations from reaching the patient, as experienced with larger, fixed-wheel systems.

Comfort. Reduce patient handling risks and enhance patient comfort. The BackSaver Fowler® feature provides the caregiver with easier head section articulation. Patients remaining on stretchers for extended periods of time due to overcrowding are at risk for developing pressure ulcers. Our innovative, five-inch air mattress—AccuMax Quantum™ Stretcher Pad—is the only non-powered, dynamic pressure redistribution stretcher mattress with Controlled Release™ Technology.

Quality. Our reinforced frames are made of heavy-gauge, automotive-grade steel for long-lasting quality. Hill-Rom is committed to helping caregivers manage an increasingly obese population. Focused on patient and caregiver safety, all Hill-Rom stretchers are designed and tested to safely support a weight of 700 lbs.

Hill-Rom

Enhancing Outcomes for Patients and Their Caregivers.™

Hill-Rom®

Procedural Stretcher

| Features | Standard Features | Optional Features |
|-----------------------------------------------------------------|-------------------|-------------------|
| 700 lb weight capacity | ● | |
| 26" width | ● | |
| 30" width | | ● |
| BackSaver Fowler™ (includes manual knee gatch) | | ● |
| Auto Contour™ option (includes "BackSaver Fowler") | | ● |
| Manual knee gatch | | ● |
| Active hand brake | | ● |
| Brake/Steer pedals – head and foot | ● | |
| Brake/Steer pedals – all four sides | | ● |
| Footend hi/lo hydraulic pedals | | ● |
| Next generation Steering Plus™ system | ● | |
| Integrated O ₂ tank holder with utility tray | | ● |
| Integrated dual display scale | | ● |
| Ergonomic push handles | | ● |
| Ergonomic push handles with integrated IV pole transport holder | | ● |
| Removable IV pole - 2 stage | | ● |
| Permanent IV pole - 3 stage | | ● |
| Siderail system (tuck-away) | ● | |
| Radiolucent deck | | ● |
| 3" density foam mattress with fire barrier | ● | |
| 4" density foam mattress with fire barrier | | ● |
| 5" Comfortline multi-density foam with fire barrier | | ● |
| 3" or 4" Tempur-Pedic | | ● |
| 5" AccuMax Quantum™ Stretcher Pad (Elastimax or Hypolex) | | ● |
| Bumper colors – choice of neutral, blue, red, teal or purple | ● | |

Accessories

- Accessory IV pole
- Ankle stirrups
- Armboard
- Chart holder
- Convertible footboard
- Detachable footboard
- Foot extender
- Infusion support system
- IV transporter
- Lateral X-ray cassette holder
- Liquid oxygen tank holder
- Oxygen tank holder
- Paper roll dispenser
- Patient tray
- Restraint straps
- Siderail covers
- Stretcher pillow
- Transport straps
- Utility tray

Please confirm accessory availability for model of stretcher ordered.

Technical Specifications

Hill-Rom® Procedural Stretcher


| | |
|----------------------------------------------|-------------------------------------|
| Maximum weight limit..... | 700 lbs (317.5 kg) |
| Low position (floor to top of surface) | 20.7" (52.6 cm) |
| with scale..... | 20.5" (52.1 cm) |
| High position | |
| floor to top of surface..... | 34.25" (87 cm) |
| with scale..... | 34" (86.4 cm) |
| Overall length | 83" (210.8 cm) |
| Overall width (siderails up) | 32" (81.3 cm) or 36" (91.4 cm) |
| Overall width (siderails stored) | 26.5" (67.3 cm) or 30.25" (76.8 cm) |
| Siderail length:..... | 47" (119.4 cm) |
| Siderail height above mattress/surface..... | 14" (35.6 cm) |
| Mattress size | 26" (66 cm) x 75" (190.5 cm) |
| | or 30" (76.2 cm) x 75" (190.5 cm) |
| Maximum head elevation | 90° |
| Maximum Trend/Reverse Trend..... | 18° |
| Floor to base clearance | 3.5" (8.9 cm) |
| Caster size..... | 8" (20.3 cm) |



AccuMax Quantum™ VPC Mattress



Tempur-Pedic® Medical Mattress

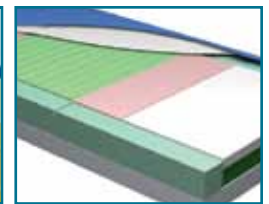

5" Comfortline® Foam Mattress



Active Hand Brake



Integrated Utility Tray



Ergonomic Push Handles



Built-in Scale



BackSaver Fowler™



O₂ Tank Holder

Hill-Rom reserves the right to make changes without notice in design, specifications and models. The only warranty Hill-Rom makes is the express written warranty extended on the sale or rental of its products.

© 2010 Hill-Rom Services, Inc. ALL RIGHTS RESERVED.

Hill-Rom

Enhancing Outcomes for Patients and Their Caregivers.™

USA 800-445-3730
Canada 800-267-2337

www.hill-rom.com